



# 1246 - High Resolution Laser spectroscopy

▷ Diode laser characteristics:

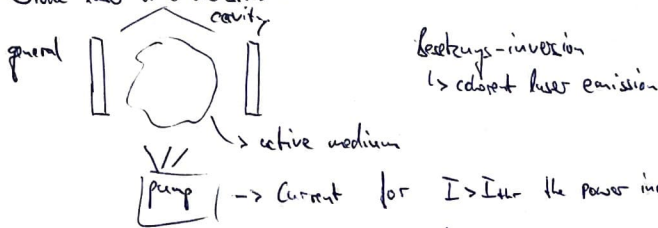


simple diode laser: p-n junction in forward bias: (homojunction)

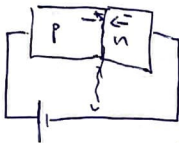

better: double heterojunction → smaller active region → high divergence and also astigmatism

LNTR configuration → smaller linewidth

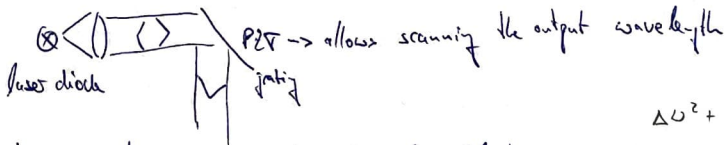


Tuning of output frequency:

- Laser temperature (slow)
- Laser current (fast)
- grating angle (only in single mode)

$$\Delta \omega^2 + \Delta \omega^2 \cdot \frac{\kappa}{P_{sat}} = 0$$

$$\Rightarrow \kappa = -P_{sat}$$

▷ Laser power and intensity:

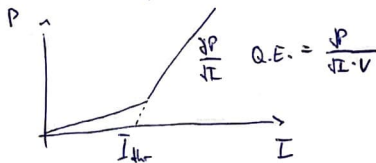

$$\sqrt{1} = e^{-\chi \frac{1}{1+s}} \cdot L$$

$$\approx 1 - \alpha \frac{1}{1+s}$$

$$R_{ing} \propto R = 1 - \sqrt{1 - \frac{\alpha}{1+s}}$$

▷ Fabry - Perot interferometer

periodic transmission:  $T_{FPI} = \frac{1}{1 + F \sin^2 \frac{\phi}{2}}$


$$\frac{L}{\lambda} = n k_0 L = \frac{n \omega L}{c}$$

$$F = \frac{4R}{(1-R)^2}$$

free spectral range:  $\Delta \omega_{FSR} = \frac{c}{2nL}$

FWHM =  $\frac{\Delta \omega_{FSR}}{F}$   $F = \frac{\pi R}{1-R}$  considering also transv. modes:  $\Delta \omega = \frac{c}{2nL} = \frac{c}{4nL}$

▷ Properties of rubidium Notation:  $n^{2s+1}L_J$

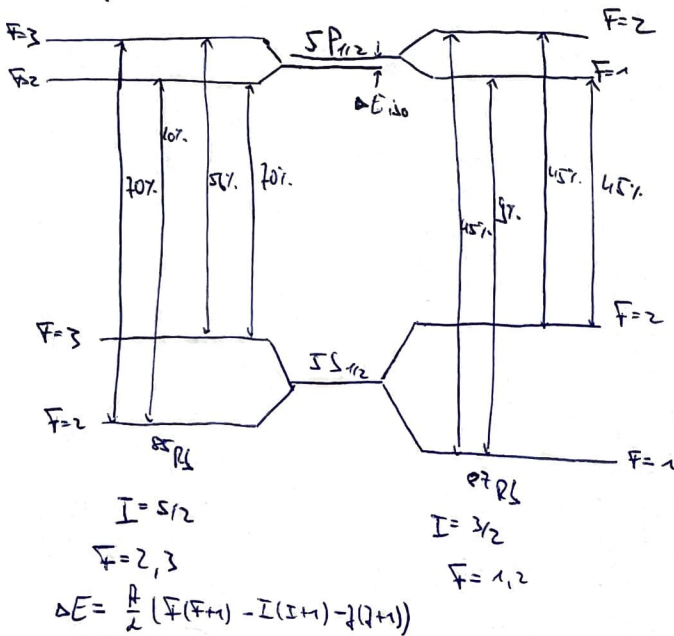

▷ linear spectroscopy

- **LORENTZ-LEER**:  $\Gamma(\omega) = e^{-\chi(\omega)L}$   $\chi(\omega) = \chi_0 \frac{(\Delta\omega/2)^2}{(\omega - \omega_0)^2 + (\Delta\omega/2)^2}$  (LORENTZ)  
↳ natural width  $\tau = \frac{1}{\Gamma} = \frac{1}{\Delta\nu_{\text{LWR}}}$
- **DOPPLER-width**:  $\chi_0(\omega) = \chi_0(\omega_0) \exp\left[-\frac{(\omega - \omega_0)^2}{2\Delta\omega_D^2}\right]$   
↳ obtained by convoluting  $\chi(\omega)$  with velocity distribution

▷ non-linear spectroscopy

- **saturation spectroscopy**: significant fraction of atoms is in excited state

$$\chi(\omega) \rightarrow \tilde{\chi}(\tilde{\omega}) \quad \Delta\tilde{\omega} = \omega \sqrt{1+s}$$

absorbed power:  $\Delta P \propto \frac{s}{1+s} \rightarrow \infty$  for  $s \rightarrow \infty$  "optical power broadening"

pump beam: turn LORENTZ-like into absorption spectrum

probe beam: measure the absorption on resonance  $\rightarrow$  LOKS-Dip

(cross-over resonances  $\rightarrow$  as name suggests)

$$E = h \cdot \omega = h \cdot \frac{c}{\lambda}$$

$$\omega = \frac{c}{\lambda}$$

$$\lambda = \frac{c}{\omega}$$

